



NRF2/ARE Secreted Luciferase Reporter Cell Line

(For Research Use Only)

HEK293 Catalog Number: SL-4042

Introduction

NRF2 plays a crucial role in cellular antioxidant defense, making it a therapeutic target for neurodegenerative diseases and cancer. Under normal conditions, NRF2 localizes in the cytosol and is rapidly degraded by the proteasome. Under oxidative stress, NRF2 is stabilized and translocates to the nucleus where it binds to a DNA promoter and initiates gene expression. In the nucleus, NRF2 forms a heterodimer with a small Maf protein and binds to the Antioxidant Response Element in the upstream promoter region of many antioxidative genes and initiates their transcription.

Product description

Signosis has developed NRF2/ARE Secreted luciferase reporter stable cell line by co-transfecting NRF2/ARE Gaussia luciferase reporter vector and hygromycin expression vector. The hygromycin-resistant clones were subsequently screened for TBHQ-induced luciferase activity. The clones with the highest fold induction were selected and expanded to produce those stable cell lines.

Materials provided


One vial of 2×10^6 cells, at passage 4, in Freezing Media. **IMPORTANT:** store the frozen cells in liquid nitrogen until you are ready to thaw and propagate them.

Handling cells upon arrival



It is strongly recommended that you propagate the cells by following instructions as soon as possible upon arrival.**

IMPORTANT: It is imperative that an adequate number of frozen stocks be made from early passages as cells may undergo genotypic changes. Possible genetic instability in transfected cells may result in a decreased responsiveness over time in normal cell culture conditions.

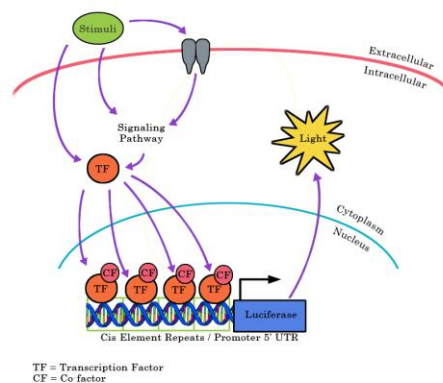

TF = Transcription Factor
CF = Co factor

Required Cell Culture Media

- **Complete Growth Media**
In 450mL of DMEM, add 50mL FBS (10% final) and 5mL Penicillin/Streptomycin (1% final). Note: HepG2 cells grow better in DMEM with 5-8% FBS.
- **2x Freezing Media**
Add 10% DMSO (final) to Complete Growth Media and sterile filter. Make fresh each time.

Materials required but not provided (Can be substituted with comparable third-party products):

Materials	Product number
Dulbecco's Modified Eagles Medium (DMEM)	Cytiva SH30243.FS
Fetal Bovine Serum (FBS)	Cytiva SH30910.03
Penicillin/Streptomycin	Cytiva P/N SV30010
Trypsin	Cytiva SH30042.02
Phosphate-buffered saline (PBS)	Cellgro P/N 21-040-CV
DMSO	Sigma P/N D8418
96-well white plate	Greiner Bio-One P/N 655098
Gaussia Luciferase substrate	Signosis P/N GLUC010
Hygromycin B (optional)	SCBT P/N sc-29067

Initial Culture Procedure

1. Quickly thaw cells in a 37°C water bath with careful agitation. Remove from the bath as soon as the vial is thawed.
2. Transfer cells to a 100mm² dish (or T-25cm² flask) containing 10ml of Complete Growth Media.
3. Gently rock the flask to ensure the cells are mixed well in the media. DO NOT PIPET.
4. Place the flask with cells in a humidified incubator at 37°C with 5% CO₂.
5. After cells adhere (wait at least 8 hours to overnight), replace media with fresh Complete Growth Media.

Subculture Procedure

1. After Cells have recovered and growing well subculture/passage cells when the density reaches 90-100% confluency, maintain and subculture the cells in Complete Growth Media.
Note: During the time that cells are not used for the experiment ideally, they can be maintained in Complete Growth Media with 50-100µg/ml of Hygromycin B.
2. Carefully remove the culture media from cells by aspiration.
3. Rinse cells with PBS, being careful to not dislodge attached cells. Then remove PBS by aspiration.
4. Add 1-2 mL trypsin/Tris-EDTA solution.
5. Incubate with trypsin for 2-5 minutes (or until detached). Confirm detachment by observation under the microscope.
6. Add 5-10ml of pre-warmed Complete Growth Media and gently pipet up and down to break the clumps.
7. Passage cells in 1:3 to 1:5 ratio when they reach 90% confluency.

NOTE: Stable cell lines may exhibit a slower proliferation rate compared to parental cells. Do not seed cells at suboptimal density as this may hinder cell growth and division.

Preparing frozen stocks

This procedure is designed for 100mm² dish or T-75cm² flasks. Scale volumes accordingly to other vessels.

1. When cells reach 90-100% confluency, freeze down cells.
2. Detach cells according to "Subculture Procedure."
3. Transfer cells to a 15ml conical centrifuge tube and centrifuge at 250 x g (or 2,000 RPM) for 5 minutes to collect the cells into a pellet.
4. Carefully aspirate the media.
5. Add 1ml of **Freezing Media** and gently resuspend by pipetting up and down.

6. Transfer 1mL of cells into a cryogenic vial.
7. Place the cryogenic vial in a freezing container (Nalgene # 5100-0001) and store it at -80°C freezer overnight.
8. Transfer cells to liquid nitrogen for long-term storage.

Assay procedure

NOTE: Gaussia Luciferase works with coelenterazine based substrates, such as our Gaussia Luciferase Substrate (cat# GLUC010). This cell line is not intended to be used with other kinds of substrates.

The following procedure should be followed as a guideline. You will need to optimize the assay conditions based on your experimental setup.

1. The day before performing the assay, trypsinize the cells and seed each well of a **white clear-bottom** 96 well plate with 1-3 x 10⁴ cells in 100µl medium.
2. Incubate the plate in a humidified incubator at 37°C with 5% CO₂ overnight.
3. Add inducing reagent directly to each well and incubate for an appropriate time to produce maximal induction.
4. Remove 10-20 µl of media from treated wells after treatment is over and plate them to a white or black opaque 96-well plate.
Note: You can collect media at different time points to monitor luciferase expression without sacrificing the cells.
5. Add 50µl of Gaussia luciferase substrate to each well and gently pipette up and down.
6. Immediately read the plate in a luminometer.

Note: We recommend a luminometer with a sensitivity of at least 3x10⁻²¹ moles luciferase.

For Data, visit

[http://www.signosisinc.com/data/Luciferase Reporter Stable Cell Lines](http://www.signosisinc.com/data/Luciferase_Reporter_Stable_Cell_Lines)

Signosis Luciferase Reporter Stable Cell Lines

For a complete list of cell lines please visit our website at <http://www.signosisinc.com/category/cell-based-assays>

Transcription Factor	Pathway	Cell Line	Cat #
NFkB	NFkB	Hela; human cervical cancer	SL-0001
NFkB	NFkB	NIH/3T3; mouse fibroblast	SL-0006
NFkB	NFkB	HEK293; human embryonic kidney	SL-0012
NFkB	NFkB	MCF-7; human breast cancer	SL-0013
NFkB	NFkB	A549; human lung cancer	SL-0014
NFkB	NFkB	HepG2; human liver cancer	SL-0017
NFkB	NFkB	MEF; murine embryonic fibroblast	SL-0033
NFAT	Calcium Signaling	Jurkat; human T lymphocytes	SL-0032
NFAT	Calcium Signaling	Hela; human cervical cancer	SL-0018
p53	p53	Hela; human cervical cancer	SL-0011
p53	p53	RKO; human colon cancer	SL-0007
SMAD	TGFbeta	HepG2; human liver cancer	SL-0016
SMAD	TGFbeta	NIH/3T3; mouse fibroblast	SL-0030
NRF2	Antioxidant Response	MCF7; human breast cancer	SL-0010
STAT1	JAK-STAT	Hela; human cervical cancer	SL-0004
STAT3	JAK-STAT	Hela; human cervical cancer	SL-0003
HIF	Hypoxia Response	NIH/3T3; mouse fibroblast	SL-0005
HIF	Hypoxia Response	Hela; human cervical cancer	SL-0023
HIF	Hypoxia Response	Neuro2a; mouse neuroblastoma	SL-0027
ER	Estrogen Receptor Signaling	T47D; human breast cancer	SL-0002
AR	Androgen Receptor Signaling	MDA-MB-453; human breast cancer	SL-0008
GR	Glucocorticoid Receptor Signaling	MDA-MB-453; human breast cancer	SL-0009
GR	Glucocorticoid Receptor Signaling	Hela; human cervical cancer	SL-0021
AP-1	JNK, ERK, MAPK Signaling	Hela; human cervical cancer	SL-0019
CREB	cAMP, PKA, CaMK Signaling	HEK293; human embryonic kidney	SL-0020
CREB	cAMP, PKA, CaMK Signaling	NIH/3T3; mouse fibroblast	SL-0031
CHOP	Unfolded Protein Response, ER stress	Mia-Paca2; human pancreatic cancer	SL-0025
TCF/LEF	Wnt/b-catenin	HEK293; human embryonic kidney	SL-0015
TCF/LEF	Wnt/b-catenin	Hela; human cervical cancer	SL-0022
TCF/LEF	Wnt/b-catenin	CHO-K1; Chinese Hamster Ovary	SL-0028
ELK	MAPK Signaling	HEK293; human embryonic kidney	SL-0040
ELK	MAPK Signaling	Hela; human cervical cancer	SL-0041
IRF	Immune Response Pathway	HEK293; human embryonic kidney	SL-0035

** Signosis products are warranted for 30 days from the date of shipment, and this warranty is valid only if the product is stored and handled according to this user manual or product information sheet and shipped directly by Signosis.